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## **SAFETY INSTRUCTIONS**

Read instructions in manual carefully before use.

This appliance may be used by children of more than 8 years of age, persons with physical, sensory or mental capabilities or inexperienced persons only if they have been previously trained in using it in a safe manner, and only if they have been informed of the dangers the product can cause.

Children may not play with the appliance. Children may not clean or carry out any maintenance on the appliance without adult supervision.

Do not use in environments with over 85% humidity.

Keep product in a clean, dry place.

Do not use accessories other than those in the instructions.

Do not handle or try to open the product (always get technical assistance).

The guarantee will be invalidated.

Do not expose product to direct sunlight or immerse in water.

Do not use in close proximity to mobile phones, televisions or microwave ovens.

Do not use batteries other than those in the instructions and do not attempt to recharge them. Remove the batteries when the product is not in use for a lengthy period.

The product is not advised for people who suffer from epilepsy or heart conditions; people who use cardiac simulators.

Do not use on wet or damp hair.

After use of comb, it is advisable to use a repellent treatment (shampoo, cream) to prevent further infestation from lice (for children in particular).


Disinfect the comb thoroughly before use by any other person or family member. To disinfect, use a damp cloth soaked in a 70% alcohol solution. Before re-use, make sure the comb is completely dry after cleaning.

Do not use if teeth show any sign of defect.

#### **USING THE PRODUCT**

The comb has been designed to fight human parasites (head lice) in a simple, fast and safe manner, with no damage to the hair or scalp. The comb is totally safe and innocuous even for very young children. When the teeth come into contact with lice, a slight electric tension is generated which is imperceptible to the human body and kills the lice. The lice are removed from the hair and scalp by the steel teeth of the comb. During use, no substance is released from the product. It is therefore possible to repeat the operation at any time. A periodic oscillatory wave is generated by the comb when it

comes into contact with the lice. A very slight electric tension is generated which instantly kills the lice. The tension generated by the included batteries is extremely slight, absolutely imperceptible to the human body and causes no harm whatsoever. Open the battery compartment and insert 1 AAA battery (making sure insertion correctly follows polarity).

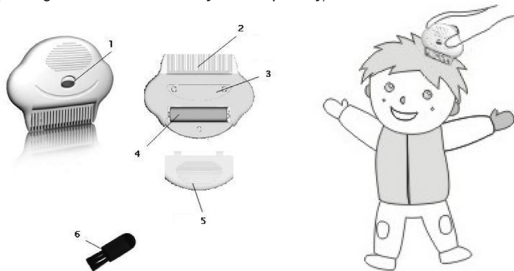- 1- On/Off button
- 2- Teeth of comb
- 3- Label
- 4- Battery compartment
- 5- Battery cover
- 6- Brush for cleaning

Press the On/Off button. A sound will indicate the comb is ready for use.

Hold the comb at a 45° angle to the scalp. The LED on the comb should be in a frontal position to the scalp. Lice usually nest in the warmer areas of the head (generally near the neck).

Move the comb slowly through the hair, preferably from zone to zone, to achieve best results. The movement of the comb should go from the back to the front of the head and then vice versa, making sure a 45° angle is maintained.

The sound will be interrupted each time a louse is picked out and killed. If the sound does not start up again in 5 seconds, use the small cleaning brush provided. Two of the teeth may come into contact during use and this may cause a sudden interruption of the sound.

Keep combing carefully and in depth, until all the lice have been killed and removed. The comb may be used at any time of the day or night.

After each use, it is advisable to use a repellent treatment (shampoo, cream) to prevent further infestation from lice (for children in particular).

During periods of more consistent lice infestation, the comb should be used frequently and preferably before going to bed.

## **TECHNICAL INFORMATION**

Frequency : 5.1 kHz

Frequency of sound emitted : 6,9 kHz

Sound indicator

Power : 1 x AAA (Test battery is included). Insert the batteries in the compartment with the poles correctly aligned. Replace dead batteries with new ones.

**The supplied batteries are for trial use only and they may have a shorter life.**

This item is made in conformity with all the applicable European directives

#### **DISPOSAL**

The device (including its removable parts and accessories) must not be disposed of together with municipal waste at the end of its life, but in compliance with European Directive. Since it must be handled separately from household waste, either carry it to a separately collected waste disposal centre for electrical and electronic appliances or give it back to the retailer on purchasing a new device with the same purpose. Any infringement will be severely prosecuted. The batteries used in this device must be disposed of in the special bins at the end of their life.

*Specifications and designs are based on the latest information available at the time of printing and subject to change without notice.*

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**JOYCARE S.p.A. socio unico**